

## Anti-Cholera toxin antibody, rabbit serum

64-007 100 ul

**Cholera toxin**, a main enterotoxin, interacts with G proteins and increases cyclic AMP in the intestinal lining to open ion channels. As ions flow into the intestinal lumen (lining), body fluids (mostly water) flows out of the body due to osmosis leading to massive diarrhea as the fluid is expelled from the body.

**Cholera toxin** is a complex consisting of one molecule of A subunit (27.2 kD) and 5 molecules of B subunits (11.6 kD). After secretion, A subunit is proteolytically processed into A1 (22 kD) and A2 (5 kD) subunits which are held together by a disulfide bond. The toxin adsorbs to GM1 ganglioside on the surface of target cells by the B subunit and the A subunit is dissociated from the B subunit during penetration. The A subunit constitutively activates adenylyl cyclase activity of  $\alpha$  subunit of Gs (a kind of GTP-binding protein).

### Applications:

1. Western blotting (dilution: 1/2,000)
2. Immunoprecipitation
3. ELISA.

Other applications have not been tested.

**Immunogen:** Cholera toxin and the toxoid purified from culture medium of *Vibrio cholerae* 569B strain

**Form:** Antiserum added with 0.05% sodium azide

**Storage:** Shipped at 4°C and stored at -20°C

**Data Link:** UniProt KB [Cholera toxin](#)

### References:

1. Hirst TR and D'Souza JM In *The Comprehensive Sourcebook of Bacterial Protein Toxins* Alouf J and Popoff M ed. 3<sup>rd</sup> edn. p. 270-290 Academic Press (2006)
2. Finkelstein RA and LoSpalluto JJ "Pathogenesis of experimental cholera. Preparation and isolation of cholera toxin and cholera toxinoid." *J. Exp Med* **130**: 185-202 (1969) PMID: [4978880](#)

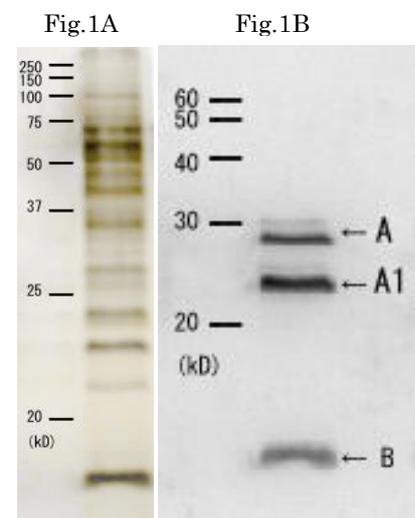


Fig.1 A SDS-PAGE analysis of culture medium of *Vibrio cholerae*. Culture medium of *Vibrio cholerae*, 569B strain was subjected to electrophoresis under reducing condition followed by silver-staining.

Fig.1 B Western blotting of culture medium of Cholera toxin.

Culture medium of *Vibrio cholerae* 569B strain was subjected to electrophoresis under reducing condition followed by Western blotting using this antibody (1/2000 dilution). A, A1 and B indicate the subunits A, A1 and B, respectively. A2 subunit (5 kD) is too small to be seen by this analysis.

**Related Products:** [#01-511](#) Cholera toxin, [#01-521](#) Cholera toxin A subunit, [#01-525](#) Cholera toxin B subunit



COSMO BIO CO., LTD.

TOKYO JAPAN

TEL : (81)3-5632-9617

FAX : (81)3-5632-9618

e-mail : [export@cosmobio.co.jp](mailto:export@cosmobio.co.jp)



[www.cosmobio.com](http://www.cosmobio.com)