



MONOCLONAL ANTIBODY

For research use only. Not for clinical diagnosis.

Catalog No. CE-036A

Anti- Histone H3 T32ph

BACKGROUND


Post-translation modifications of histones modulate the accessibility and transcriptional competence of specific chromatin regions within the eukaryotic genome. Phosphorylation of histone H3 is unique in the sense that it associates on one hand with open chromatin during gene activation and marks on the other hand highly condensed chromatin during mitosis.

Product type	Primary antibodies
Immunogen	Synthetic peptide corresponding to N-terminus region Thr32ph (aa 21-39) of human Histone H3, ATKAARKSAPS(phT)GGVKKPH
Rased in	Rat
Myeloma	SP2
Clone number	6G8B7
Isotype	IgG2a, κ
Host	-
Source	Culture supernatant
Purification	Ion-exchange chromatography
Form	Liquid
Storage buffer	PBS containing 50% Glycerol, 0.05% NaN_3 as a preservative
Concentration	1 mg / ml
Volume	100 μl
Label	Unlabeled
Specificity	Histone H3 T32ph Epitope : phosphorylated Thr32 of Histone H3
Cross reactivity	Human, Monkey, Mouse, Rat, Hamster Other species have not been tested.
Storage	Store below -20°C (below -70°C for prolonged storage) Aliquot to avoid cycles of freeze/thaw.
Other	Data Link : UniProtKB/Swiss-Prot P68431 * recommended positive controls is mammalian cell

Application notes	• Western blotting: 1/500 (Fig.2)
Recommended dilutions	• Immunocytochemistry: 1/500 (Fig.3)

Other applications have not been tested.
Optimal dilutions/concentrations should be determined by the end user.

References	1) Strahl and Allis, (2000) Nature 403, 41-45.
-------------------	--

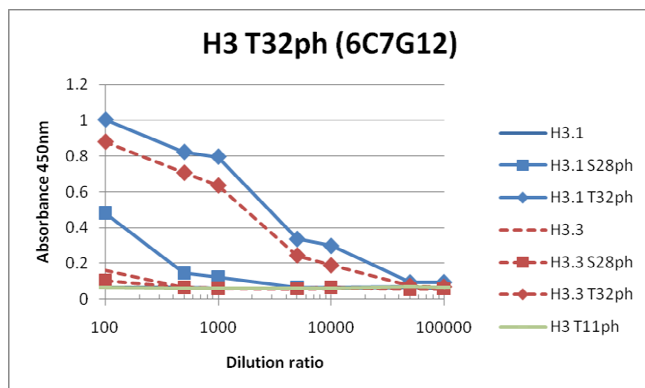

H3.1	21	ATKAARKS	APAT	GGVKKPH	39
H3.1 S28ph	21	ATKAARKS	phAPAT	GGVKKPH	39
H3.1 T32ph	21	ATKAARKS	APATph	GGVKKPH	39
H3.3	21	ATKAARKS	APST	GGVKKPH	39
H3.3 S28ph	21	ATKAARKS	phAPST	GGVKKPH	39
H3.3 T32ph	21	ATKAARKS	APSTph	GGVKKPH	39
H3 T11ph	1	ARTKQTARKST	phGGKAPRKQC	19	

Fig.1 The composition of Histone H3 peptides and the reactivity of Histone H3 T32ph antibody, 6C7G12.

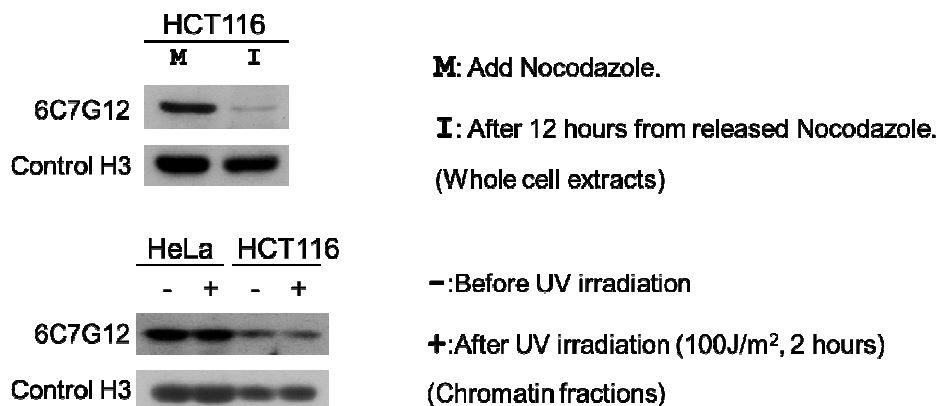


Fig.2 Western blot analysis of the treated-cell extracts using Histone H3 T32ph antibody, 6C7G12.

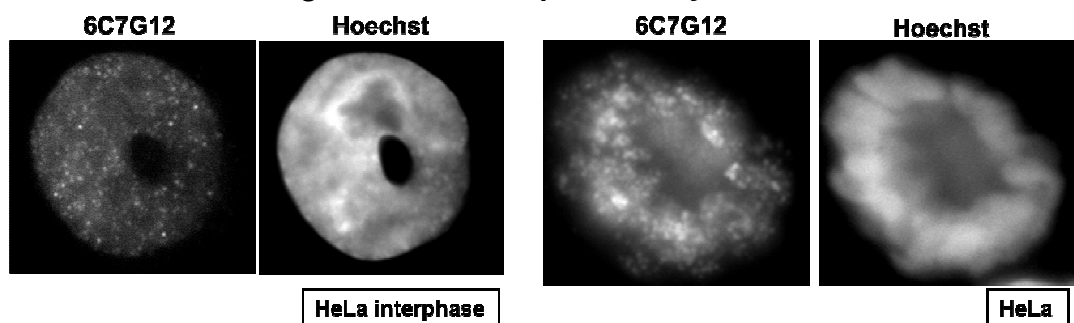


Fig.3 Immunocytochemical analysis of eLa Cell using Histone H3 T32ph antibody, 6C7G12. (left : interphase, right : metaphase).

RELATED PRODUCTS:

- Anti Histone H3 S10ph Monoclonal Antibody [Cat# CE-034A]
- Anti Histone H3 T11ph Monoclonal Antibody [Cat# CE-035A]
- Anti Histone H3 T32ph Monoclonal Antibody [Cat# CE-036A]
- Anti Histone H3 K9Ac Monoclonal Antibody [Cat# CE-037A]
- Anti Histone H3.1 Monoclonal Antibody [Cat# CE-039A]
- Anti Histone H3.1 Monoclonal Antibody [Cat# CE-039B]
- Anti Histone H3.3 Monoclonal Antibody [Cat# CE-040A]
- Anti Histone H3.3 Monoclonal Antibody [Cat# CE-040B]

For research use only. Not for clinical diagnosis.



COSMO BIO CO., LTD.

Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

http://www.cosmobio.co.jp/index_e.asp

E-mail: export@cosmobio.co.jp

Phone : +81-3-5632-9617

FAX : +81-3-5632-9618

New

World's First !!!

Histone Variant Monoclonal Antibodies

Anti Histone H3.1/H3.2 [Clone: 6G3C7]

Anti Histone H3.1/H3.2 [Clone: 1D4F2]

Anti Histone H3.3 [Clone: 6C4A3]

Anti Histone H3.3 [Clone: 1E4A3]


Nucleosomes are composed of four different histone proteins designated H2A, H2B, H3, and H4. In humans, five variants of histone H3 are reported: H3.1, H3.2, H3.3, H3t, and CENP-A. The two major Histone H3 variants, H3.1 and H3.3, are the main variants displaying distinct genomic localization patterns in eukaryotes. Deposition of Histone H3.1 is associated with DNA synthesis during DNA replication and possibly DNA repair, while Histone H3.3 is incorporated independently of DNA synthesis and is the predominant form of H3 found in non-dividing cells. Hence, these new Histone H3 variant monoclonal antibodies

offer great utility for dissecting the functional significance of these H3 variants and the molecular mechanisms associated with their deposition.


Recently, it was shown that a genomic gene cluster regulating skeletal myogenesis is marked by H3.3 protein prior to cellular muscle formation and that H3.3 marking of this region enables myogenic gene activation (Ref. 2). These results suggest that monitoring H3.3 marking at specific loci may be useful in the prediction of cell fate. These H3.3 monoclonal antibodies are expected to be useful probes in the field of regenerative medicine.

Antibody specificity by competition peptide ELISA


Histone H3.1/H3.2 MAb (6G3C7)

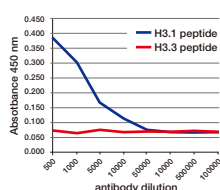


Histone H3.1/H3.2 MAb (1D4F2)

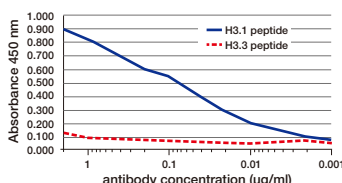

Histone H3.3 MAb (6C4A3)

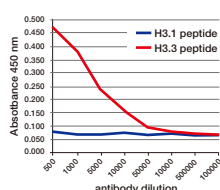

Histone H3.3 MAb (1E4A3)

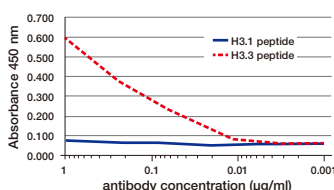


H3.1 peptide 79 KTDLRFS~~SA~~VMALQEASEA 97
H3.3 peptide 79 KTDLRFS~~AA~~IGALQEASEA 97

H3.1 peptide 21 ATKAARKSAP~~AT~~GGVKKPH 39
H3.3 peptide 21 ATKAARKSAP~~ST~~GGVKKPH 39


Fluorescence immunostaining

Histone H3.1/H3.2 MAb (1D4F2)

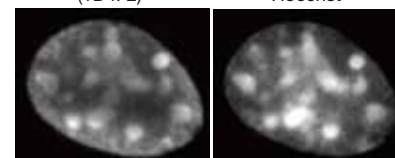

NIH3T3

Histone H3.3 MAb (1E4A3)

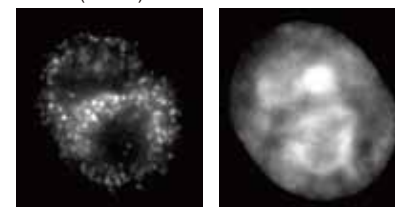

HeLa

Experimental example

These H3 variant antibodies were essential tools in a first of kind study showing that differentiation specific genes are marked for lineage specific expression by the deposition of Histone H3.3 at the onset of differentiation signaling (Ref. 2).

Reference


- 1) Hake and Allis, (2006) *PNAS*, 103, 6428-6435.
- 2) Harada et. al., (2012) *EMBO J.* 36, 2994-3007.

H3.3 gene marking in skeletal muscle differentiation

MARKING



Non-MARKING

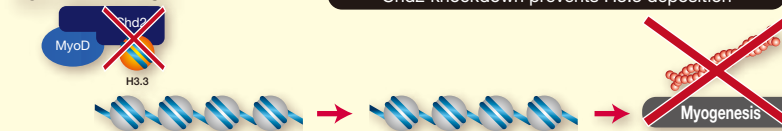

Description	Clone	Isotype	Epitope	Application	Cat. No.	Quantity
Anti Histone H3.1/H3.2	6G3C7	Rat-IgG1, λ	H3.1/H3.2 (79-94)	IP/ WB	CAC-CE-039A	100 μ L (100 μ g)
Anti Histone H3.1/H3.2	1D4F2	Mouse-IgG2b, λ	H3.1/H3.2 (21-39)	ChIP/ IP/ WB/ IC	CAC-CE-039B	50 μ L (50 μ g)
Anti Histone H3.3	6C4A3	Rat-IgG2a, κ	H3.3 (79-97)	IP/ WB	CAC-CE-040A	100 μ L (100 μ g)
Anti Histone H3.3	1E4A3	Rat-IgG2a, λ	H3.3 (21-39)	ChIP/ IP/ WB/ IC	CAC-CE-040B	50 μ L (50 μ g)




COSMO BIO CO., LTD.


Monoclonal Antibodies to Histone Modifications

Histones are the main protein components of chromatin. To facilitate nuclear packaging and control of gene expression, DNA in chromatin is wound around nucleosome particles composed primarily of the Histones H2A, H2B, H3, and H4. Histone N-terminal regions (histone tails) protrude from the nucleosome core and are subject to a variety of reversible, regulated modifications (including acetylation, phosphorylation, and methylation) influencing transcription and chromatin structure. How such modifications are regulated and how these modifications effect gene expression continues to be an area of intense interest and research. In such studies, chromatin immunoprecipitation (ChIP) is perhaps the most widely used experimental procedure. Due to the inherent variability and limited supply of polyclonal antibodies, well characterized monoclonal antibodies are preferred reagents for ChIP. The versatile set of anti-histone monoclonal antibodies offered here are therefore highly valuable reagents to your lab's epigenetic toolbox.

Histone H3 N-terminal modifications

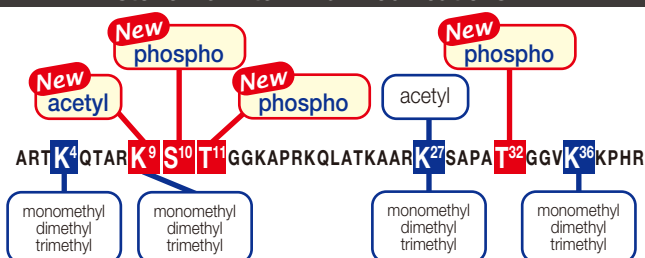

Histone H3 phospho Ser10 immunostaining

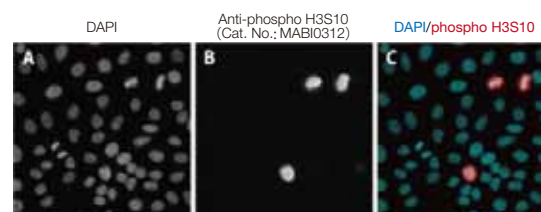

Description	Host	Residue	Modification	Clone	Application	Cat. No.	Quantity
Anti Histone H3	Mouse	-	unmodified	MAB10301	ChIP/ WB/ IC	MCA-MAB10001-100-EX	100 µL (100 µg)
Anti Monomethyl Histone H3 (Lys4)	Mouse	K4 (Lysine 4)	monomethyl	MAB10302	ChIP/ WB/ IC	MCA-MAB10002-100-EX	100 µL (100 µg)
Anti Dimethyl Histone H3 (Lys4)	Mouse		dimethyl	MAB10303	ChIP/ WB/ IC	MCA-MAB10003-100-EX	100 µL (100 µg)
Anti Trimethyl Histone H3 (Lys4)	Mouse		trimethyl	MAB10304	ChIP/ WB/ IC	MCA-MAB10004-100-EX	100 µL (100 µg)
Anti Histone H3 K9Ac New	Rat	K9 (Lysine 9)	acetyl	2G1F9	ChIP/ WB/ IC/ IHC	CAC-CE-037A	100 µL (100 µg)
Anti Acethyl Histone H3 (Lys9)	Mouse		acetyl	MAB10305	ChIP/ WB/ IC	MCA-MAB10005-100-EX	100 µL (100 µg)
Anti Monomethyl Histone H3 (Lys9)	Mouse		monomethyl	MAB10306	ChIP/ WB/ IC	MCA-MAB10006-100-EX	100 µL (100 µg)
Anti Dimethyl Histone H3 (Lys9)	Mouse		dimethyl	MAB10307	ChIP/ WB/ IC	MCA-MAB10007-100-EX	100 µL (100 µg)
Anti Trimethyl Histone H3 (Lys9)	Mouse		trimethyl	MAB10308	ChIP/ WB/ IC	MCA-MAB10008-100-EX	100 µL (100 µg)
Anti Acetyl Histone H3 (Lys9/27)	Mouse	K9/27 (Lysine 9/27)	acetyl	MAB10310	ChIP/ WB/ IC	MCA-MAB10010-100-EX	100 µL (100 µg)
Anti Acetyl Histone H3 (Lys27)	Mouse	K27 (Lysine 27)	acetyl	MAB10309	ChIP/ WB/ IC	MCA-MAB10009-100-EX	100 µL (100 µg)
Anti Monomethyl Histone H3 (Lys27)	Mouse		monomethyl	MAB10321	ChIP/ WB/ IC	MCA-MAB10321-100-EX	100 µL (100 µg)
Anti Dimethyl Histone H3 (Lys27) coming soon!	Mouse		dimethyl	MAB10322	ChIP/ WB/ IC	MCA-MAB10322-100-EX	100 µL (100 µg)
Anti Trimethyl Histone H3 (Lys27)	Mouse		trimethyl	MAB10323	ChIP/ WB/ IC	MCA-MAB10323-100-EX	100 µL (100 µg)
Anti Monomethyl Histone H3 (Lys36)	Mouse	K36 (Lysine 36)	monomethyl	MAB10331	ChIP/ WB/ IC	MCA-MAB10331-100-EX	100 µL (100 µg)
Anti Dimethyl Histone H3 (Lys36)	Mouse		dimethyl	MAB10332	ChIP/ WB/ IC	MCA-MAB10332-100-EX	100 µL (100 µg)
Anti Trimethyl Histone H3 (Lys36)	Mouse		trimethyl	MAB10333	ChIP/ WB/ IC	MCA-MAB10333-100-EX	100 µL (100 µg)
Anti Histone H3 S10ph New	Rat	S10 (Serine 10)	phospho	6G8B7	WB/ IC	CAC-CE-034A	100 µL (100 µg)
Anti phospho Histone H3 (Ser10)	Mouse		phospho	MAB10312	ChIP/ WB/ IC	MCA-MAB10012-100-EX	100 µL (100 µg)
Anti Histone H3 T11ph New	Rat	T11 (Threonine 11)	phospho	6G12C5	WB/ IC	CAC-CE-035A	100 µL (100 µg)
Anti Histone H3 T32ph New	Rat	T32 (Threonine 32)	phospho	6C7G12	WB/ IC	CAC-CE-036A	100 µL (100 µg)
Anti phospho Histone H2B (Ser14)	Mouse	S14 (Serine 14)	phospho	MAB10251	ChIP/ WB/ IC	MCA-MAB10251-100-EX	100 µL (100 µg)

Reference

- 1) Strahl and Allis, (2000) *Nature* 403, 41-45.
- 2) Shimada et. al., (2008) *Cell* 132, 221-232.
- 3) Kimura H, et. al., (2008) *Cell Struct Funct.*, 33, 61
- 4) Ohhata T, et. al., (2008) *Development.*, 135, 227
- 5) Luco RF, et. al., (2010) *Science.*, 327, 996 (2010)
- 6) Rechtsteiner A, et. al., (2010) *PLoS Genet.*, 6, e1001091
- 7) Furuhashi H, et. al., (2010) *Epigenetics Chromatin.*, 3, 15
- 8) Matsui T, et. al., (2010) *Nature.*, 464, 927

For research use only. Not for diagnostic use.



COSMO BIO Co., LTD.

TOYO EKIMAE BLDG. 2-20, TOYO 2-CHOME,
KOTO-KU, TOKYO 135-0016, JAPAN
TEL : (81)3-5632-9617
FAX : (81)3-5632-9618
e-mail : export@cosmobio.co.jp
URL : www.cosmobio.com

New

世界初! コスモ・バイオから新発売!

ヒストンH3 バリエーションモノクローナル抗体

Anti Histone H3.1/H3.2 [Clone: 6G3C7]

Anti Histone H3.1/H3.2 [Clone: 1D4F2]

Anti Histone H3.3 [Clone: 6C4A3]

Anti Histone H3.3 [Clone: 1E4A3]


ヌクレオソームは H2A、H2B、H3 および H4 の 4 種のヒストンタンパク質で構成されています。また、ヒストンにはバリエーションが存在し、特にヒト由来のヒストン H3 には、H3.1、H3.2、H3.3、H3t、および CENP-A の 5 種類のバリエーションが報告されています。主なバリエーションである H3.1、および H3.3 は真核生物においては異なるゲノム局在パターンを示しています。Histone H3.1 は転写因子および転写補因子、DNA 複製因子、DNA 修復タンパク質への関与が、ヒストン H3.3 は DNA 合成に依存しないデポジションパスウェイにおいて、交換バリエーションとしての

機能が知られております。新製品ヒストン H3.1 抗体、ヒストン H3.3 抗体は、その機能解析に有用です。


また、最近、ゲノム上の骨格筋形成にかかわる遺伝子群は、細胞が筋肉形成される以前に H3.3 タンパク質で予めマーキングされ、このマーキングが形成されることで、細胞が筋肉組織を形成する能力を獲得することが、H3.3 抗体の使用によって明らかになりました (参考文献 2)。再生医療研究にも有用な本抗体を是非お試しください。

抗体の特異性を ELISA 法で確認


Histone H3.1/H3.2 MAb (6G3C7)

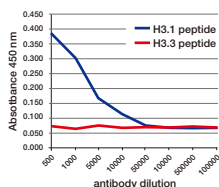


Histone H3.1/H3.2 MAb (1D4F2)

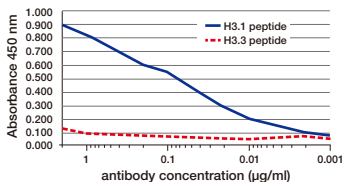

Histone H3.3 MAb (6C4A3)

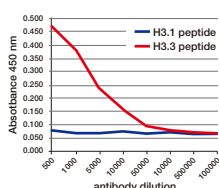

Histone H3.3 MAb (1E4A3)

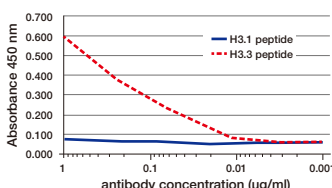



H3.1 peptide 79 KTDLRFSQSAVMALQEASEA 97
H3.3 peptide 79 KTDLRFSQSAIGALQEASEA 97

H3.1 peptide 21 ATKAAARKSAPATGGVKKPH 39
H3.3 peptide 21 ATKAAARKSAPSTGGVKKPH 39

細胞の蛍光染色

Histone H3.1/H3.2 MAb (1D4F2)



Hoechst

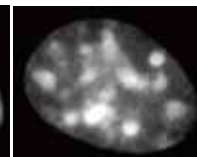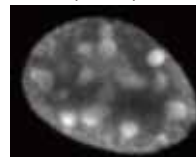

NIH3T3 細胞

Histone H3.3 MAb (1E4A3)

Hoechst

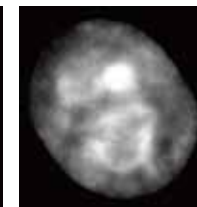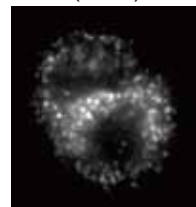

HeLa 細胞


H3.3 抗体での実験例

ES 細胞、幹細胞など特定の組織に変化する能力 (分化能) を持った細胞は、あらかじめ H3.3 と呼ばれる特殊なヒストン (ヒストンバリエーション) により、必要な遺伝子がマーキングされることが明らかになりました。骨格筋形成時に、転写因子 MyoD が骨格筋遺伝子に選択的に結合し、MyoD に結合したシャペロンである Chd2 タンパク質が H3.3 をゲノムに組み込むことでマーキングしていることが発見されました。一方で、マーキング機構を失った細胞では、骨格筋の形成がなされなかったことから、骨格筋形成を事前に予測できる細胞内のメカニズムを明らかにした世界初の知見となりました (参考文献 2)。

九州大学大学院医学研究院 准教授 大川 恭行 先生

遺伝子「マーキング」により骨格筋となる

マーキングがある場合

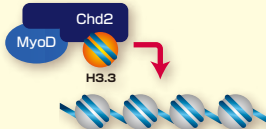

Chd2-MyoD は遺伝子を H3.3 でマーキングする



骨格筋形成

マーキングがない場合



骨格筋形成なし

参考文献 1) Hake and Allis, (2006) *PNAS*, 103, 6428-6435.
2) Harada et al., (2012) *EMBO J*, 31, 2994-3007.

コスモ・バイオ株式会社 メーカー略号: CAC

品名	クローン	アイソタイプ	エピトープ	適用	品番	包装	希望販売価格
Anti Histone H3.1/H3.2	6G3C7	Rat-IgG1, λ	H3.1/H3.2 (79-94)	WB	CE-039A	100 μ L (100 μ g)	¥ 58,000
Anti Histone H3.1/H3.2	1D4F2	Mouse-IgG2b, λ	H3.1/H3.2 (21-39)	ChIP/ WB/ IC	CE-039B	50 μ L (50 μ g)	¥ 58,000
Anti Histone H3.3	6C4A3	Rat-IgG2a, κ	H3.3 (79-97)	WB	CE-040A	100 μ L (100 μ g)	¥ 58,000
Anti Histone H3.3	1E4A3	Rat-IgG2a, λ	H3.3 (21-39)	ChIP/ WB/ IC	CE-040B	50 μ L (50 μ g)	¥ 58,000




人と科学のステキな未来へ

コスモ・バイオ株式会社

売れています

ChIP、免疫染色、WB グレード 修飾ヒストン H3/H2B モノクローナル抗体

Histone H3 N-terminal modifications

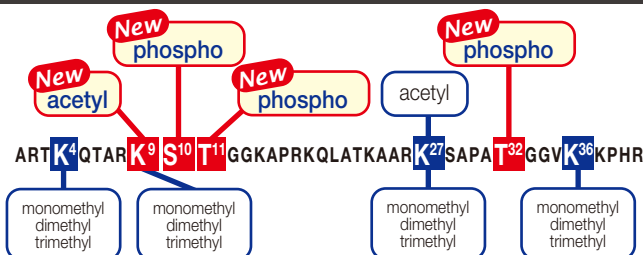

ヌクレオソームのヒストンタンパク質（H2A, H2B, H3, H4）は、クロマチンの主要な構成要素になっています。ヒストンの N 末端はヒストンテールとよばれヌクレオソームコアから少し離れて存在しています。この N 末端の部分で様々な修飾を受け細胞のシグナル刺激に応答し、アセチル化、リン酸化、メチル化など様々な修飾を受け遺伝子発現に影響を与えます。エピジェネティクス研究での ChIP アッセイに不可欠な修飾ヒストン抗体を、是非コスモラインアップからご活用ください。

M で始まる品番：株式会社モノクローナル抗体研究所
C で始まる品番：コスモ・バイオ株式会社

メーカー略号：MCA
メーカー略号：CAC

品名	免疫動物	残基	修飾	クローン	適用	品番	包装	希望販売価格
Anti Histone H3	Mouse	-	unmodified	MAB10301	ChIP/ WB/ IC	MAB10001-20	20 µL (20 µg)	¥ 18,000
						MAB10001-100	100 µL (100 µg)	¥ 58,000
Anti Monomethyl Histone H3 (Lys4)	Mouse	K4 (Lysine 4)	monomethyl	MAB10302	ChIP/ WB/ IC	MAB10002-20	20 µL (20 µg)	¥ 18,000
						MAB10002-100	100 µL (100 µg)	¥ 58,000
Anti Dimethyl Histone H3 (Lys4)	Mouse		dimethyl	MAB10303	ChIP/ WB/ IC	MAB10003-20	20 µL (20 µg)	¥ 18,000
						MAB10003-100	100 µL (100 µg)	¥ 58,000
Anti Trimethyl Histone H3 (Lys4)	Mouse		trimethyl	MAB10304	ChIP/ WB/ IC	MAB10004-20	20 µL (20 µg)	¥ 18,000
						MAB10004-100	100 µL (100 µg)	¥ 58,000
Anti Histone H3 K9Ac New	Rat	K9 (Lysine 9)	acetyl	2G1F9	ChIP/ WB/ IC/ IHC	CE-037A	100 µL (100 µg)	¥ 58,000
Anti Acethyl Histone H3(Lys9)	Mouse		acetyl	MAB10305	ChIP/ WB/ IC	MAB10005-20	20 µL (20 µg)	¥ 18,000
						MAB10005-100	100 µL (100 µg)	¥ 58,000
Anti Monomethyl Histone H3 (Lys9)	Mouse		monomethyl	MAB10306	ChIP/ WB/ IC	MAB10006-20	20 µL (20 µg)	¥ 18,000
						MAB10006-100	100 µL (100 µg)	¥ 58,000
Anti Dimethyl Histone H3 (Lys9)	Mouse		dimethyl	MAB10307	ChIP/ WB/ IC	MAB10007-20	20 µL (20 µg)	¥ 18,000
						MAB10007-100	100 µL (100 µg)	¥ 58,000
Anti Trimethyl Histone H3 (Lys9)	Mouse		trimethyl	MAB10308	ChIP/ WB/ IC	MAB10008-20	20 µL (20 µg)	¥ 18,000
						MAB10008-100	100 µL (100 µg)	¥ 58,000
Anti Acetyl Histone H3 (Lys9/27)	Mouse	K9/27 (Lysine 9/27)	acetyl	MAB10310	ChIP/ WB/ IC	MAB10010-20	20 µL (20 µg)	¥ 18,000
						MAB10010-100	100 µL (100 µg)	¥ 58,000
Anti Acetyl Histone H3 (Lys27)	Mouse	K27 (Lysine 27)	acetyl	MAB10309	ChIP/ WB/ IC	MAB10009-20	20 µL (20 µg)	¥ 18,000
						MAB10009-100	100 µL (100 µg)	¥ 58,000
Anti Monomethyl Histone H3 (Lys27)	Mouse		monomethyl	MAB10321	ChIP/ WB/ IC	MAB10321-20	20 µL (20 µg)	¥ 18,000
						MAB10321-100	100 µL (100 µg)	¥ 58,000
Anti Dimethyl Histone H3 (Lys27) 近日発売	Mouse		dimethyl	MAB10322	ChIP/ WB/ IC	MAB10322-20	20 µL (20 µg)	¥ 18,000
						MAB10322-100	100 µL (100 µg)	¥ 58,000
Anti Trimethyl Histone H3 (Lys27)	Mouse		trimethyl	MAB10323	ChIP/ WB/ IC	MAB10323-20	20 µL (20 µg)	¥ 18,000
						MAB10323-100	100 µL (100 µg)	¥ 58,000
Anti Monomethyl Histone H3 (Lys36)	Mouse	K36 (Lysine 36)	monomethyl	MAB10331	ChIP/ WB/ IC	MAB10331-20	20 µL (20 µg)	¥ 18,000
						MAB10331-100	100 µL (100 µg)	¥ 58,000
Anti Dimethyl Histone H3 (Lys36)	Mouse		dimethyl	MAB10332	ChIP/ WB/ IC	MAB10332-20	20 µL (20 µg)	¥ 18,000
						MAB10332-100	100 µL (100 µg)	¥ 58,000
Anti Trimethyl Histone H3 (Lys36)	Mouse		trimethyl	MAB10333	ChIP/ WB/ IC	MAB10333-20	20 µL (20 µg)	¥ 18,000
						MAB10333-100	100 µL (100 µg)	¥ 58,000
Anti Histone H3 S10ph New	Rat	S10 (Serine 10)	phospho	6G8B7	WB/ IC	CE-034A	100 µL (100 µg)	¥ 58,000
Anti phospho Histone H3 (Ser10)	Mouse		phospho	MAB10312	ChIP/ WB/ IC	MAB10012-20	20 µL (20 µg)	¥ 18,000
						MAB10012-100	100 µL (100 µg)	¥ 58,000
Anti Histone H3 T11ph New	Rat	T11 (Threonine 11)	phospho	6G12C5	WB/ IC	CE-035A	100 µL (100 µg)	¥ 58,000
Anti Histone H3 T32ph New	Rat	T32 (Threonine 32)	phospho	6C7G12	WB/ IC	CE-036A	100 µL (100 µg)	¥ 58,000
Anti phospho Histone H2B (Ser14)	Mouse	S14 (Serine 14)	phospho	MAB10251	ChIP/ WB/ IC	MAB10251-20	20 µL (20 µg)	¥ 18,000
						MAB10251-100	100 µL (100 µg)	¥ 58,000

お願い および 注意事項

- 希望販売価格・・・「希望販売価格」は参考であり、販売店様からの販売価格ではございません。
記載の希望販売価格は2012年11月1日現在の希望販売価格です。
予告なしに改定される場合がありますので、ご注文の際にご確認下さい。消費税は含まれておりません。
- 使用範囲・・・記載の商品は全て、「研究用試薬」です。
人や動物の医療用・臨床診断用・食品用等としては使用しないよう、十分ご注意ください。

取扱店



人と科学のステキな未来へ

コスモ・バイオ株式会社

〒135-0016 東京都江東区東陽 2-2-20 東陽駅前ビル
URL : <http://www.cosmobio.co.jp/>

● 営業部（お問い合わせ）
TEL : (03) 5632-9610 FAX : (03) 5632-9619
TEL : (03) 5632-9620

(11761)