


# EGFR (SP84)

## Rabbit Monoclonal Antibody

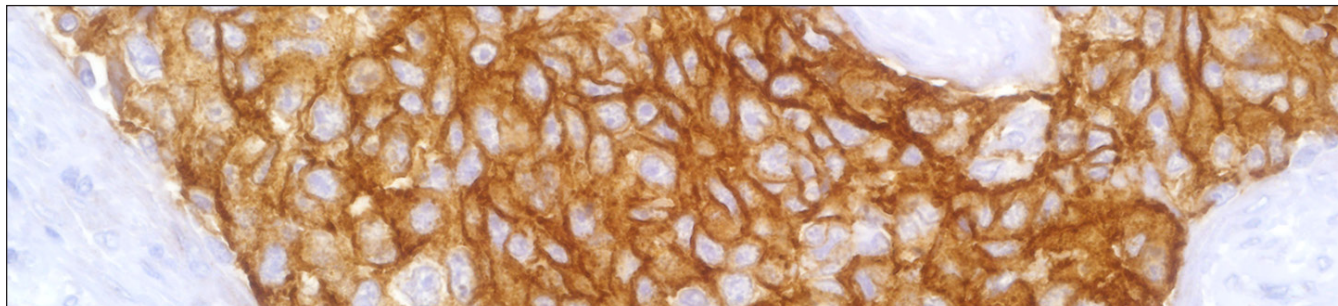

### PRODUCT AVAILABILITY

Cat. No.	Description
414R-14	0.1 ml, concentrate
414R-15	0.5 ml, concentrate
414R-16	1.0 ml, concentrate
414R-17	1.0 ml, predilute ready-to-use
414R-18	7.0 ml, predilute ready-to-use
414S	Positive control slides, 5 slides/pack

### SYMBOL DEFINITIONS

<b>P</b>	predilute	<b>E</b>	serum
<b>C</b>	concentrate	<b>DIL</b>	concentrate dilution range
<b>A</b>	ascites	<b>KEY-CODE</b>	keycode
<b>S</b>	supernatant		

### INTENDED USE

This antibody is intended for *in vitro* diagnostic (IVD) use.

The Cell Marque EGFR (SP84) antibody is intended for qualified laboratories to qualitatively identify by light microscopy the presence of associated antigens in sections of formalin-fixed, paraffin-embedded tissue sections using IHC test methods. Use of this antibody is indicated, subsequent to clinical differential diagnosis, as an aid in the identification of carcinomas within the context of an antibody panel, the patient's clinical history, and other diagnostic tests evaluated by a qualified pathologist.

### SUMMARY AND EXPLANATION

EGFR is a 170-kDa transmembrane glycoprotein encoded by the HER-1 proto-oncogene located at 7p11.2-p12.1-2, EGFR is widely expressed on the surface of epithelial cells, fibroblasts, gliocytes, keratinocytes, and

other cell types.<sup>1</sup> EGFR is overexpressed in many epithelial malignancies including carcinomas of the colorectum, stomach, esophagus, pancreas, oropharynx, adrenocortical carcinoma, non-small cell carcinoma of the lung, cutaneous and anal squamous carcinoma, and head and neck squamous carcinoma.<sup>1-5</sup> EGFR protein expression has also been a common finding in breast carcinoma, particularly in triple-negative, basal-like breast carcinomas.<sup>4</sup> Studies suggest that EGFR expression is not unique to carcinomas and may be present in malignant bone and soft tissue tumors.<sup>5</sup> Soft tissue sarcomas such as synovial sarcoma and epithelioid sarcoma show morphologic and immunophenotypic features of epithelial differentiation. Hence, EGFR overexpression is often seen in synovial sarcoma and epithelioid sarcoma. IHC analysis of 48 synovial sarcoma specimens representing primary and metastatic lesions using the anti-EGFR antibody demonstrated positive reactions in 34 of 48 cases (71%). The same study included 32 cases of malignant peripheral nerve sheath tumor, in which EGFR overexpression was found in 20 cases (62.5%).<sup>4</sup> Cascio, MJ et al. found 13 of 15 cases (87%) of epithelioid sarcoma displayed immunoreactivity of EGFR by IHC. Findings included strong, homogenous staining in the majority of cases, but absence of either gene amplification or kinase domain mutations.<sup>5</sup>

### PRINCIPLES AND PROCEDURES

The stated primary antibody (414R series) may be used as the primary antibody for immunohistochemical staining of formalin-fixed, paraffin-embedded tissue sections. In general, immunohistochemical staining in conjunction with a streptavidin-biotin detection system allows the visualization of antigens via the sequential application of a specific antibody (primary antibody) to the antigen, a secondary antibody (link antibody) to the primary antibody, an enzyme complex and a chromogenic substrate with interposed washing steps. Alternatively, a biotin-free polymer detection system may be used. The enzymatic activation of the chromogen results in a visible reaction product at the antigen site. The specimen may then be counterstained and a coverslip applied. Results are interpreted using a light microscope and aid in the differential diagnosis of pathophysiological processes, which may or may not be associated with a particular antigen.

Cell Marque prediluted products are optimally diluted for use with Cell Marque detection kits though they are commonly and successfully used with a wide variety of detection kits offered by other manufacturers.

## MATERIALS AND METHODS

See product label for lot specific information for the following:

1. Antibody immunoglobulin concentration
2. Source details

### Reagents Provided

**Prediluted** The stated primary antibody product (414R-17, 414R-18) contains ready-to-use reagent.

The predilute immunoglobulin concentration range for this product is 0.01 - 1 µg/ml.

**Concentrated** The stated primary antibody product (414R-14, 414R-15, 414R-16) contains concentrated reagent.

Both the prediluted and concentrated formats of this antibody are diluted in Tris Buffer, pH 7.3-7.7, with 1% BSA and <0.1% Sodium Azide.

The concentrate immunoglobulin concentration range for this product is 0.25 - 25 µg/ml.

The recommended working dilution range for the concentrate product is **1:25 - 1:100** and can be found on the product label.

**Isotype:** IgG

### Reconstitution, Mixing, Dilution, Titration

Prediluted antibody is ready-to-use and optimized for staining. No reconstitution, mixing, dilution, or titration is required. The concentrated antibody is optimized to be diluted to within the dilution range.

The user must validate the working dilution of the concentrated product. Differences in tissue processing and technical procedures in the laboratory may produce significant variability in results and consequently require regular use of controls. (See Quality Control Procedures section)

### Materials and Reagents Needed But Not Provided

The following reagents and materials may be required for staining but are not provided with the primary antibody:

- |  |   |
|--|---|
| 1. Positive and negative control tissue  | 9. Electric Pressure Cooker (cat. #976L) for tissue pretreatment step               |
| 2. Microscope slides, positively charged   | 10. Detection system (such as cat. #954D-20) and chromogen (such as cat. #957D-20)  |
| 3. Drying oven capable of maintaining a temperature of 58-60°C ± 5°C   | 11. Wash Solutions (cat. #935B-09)  |
| 4. Staining jars or baths  | 12. Hematoxylin (cat. #930B-05) or other counterstain                               |
| 5. Timer   | 13. Antibody diluents (such as cat. #938B-05)                                       |
| 6. Xylene or xylene substitute   | 14. Peroxide Block (cat. #925B-05) for use with HRP                                 |
| 7. Ethanol or reagent alcohol<br><i>Note: Cell Marque's one-step pretreatment, Trilogy™ (cat. #920P-06), can replace both 6 and 7 above.</i> | 15. Avidin-Biotin Block (cat. # 928B-02 for use with streptavidin-biotin detection) |
| 8. Deionized or distilled water  |   |

16. Negative Control Reagent (cat. #932B-02 for mouse; cat. #933B-02 for rabbit)

17. Mounting medium (cat. #931B-03)

18. Cover glass
19. Light microscope (40-400x)

### Storage and Handling

Store at 2-8°C. Do not freeze.

To ensure proper reagent delivery and stability of the antibody after every run, the cap must be replaced and the bottle must be immediately placed in the refrigerator in an upright position.

Every antibody reagent is expiration dated. When properly stored, the reagent is stable to the date indicated on the label. Do not use reagent beyond the expiration date for the prescribed storage method.

There are no definitive signs to indicate instability of this product; therefore, positive and negative controls should be run simultaneously with unknown specimens. Contact Cell Marque customer service if there is a suspected indication of reagent instability.

### Specimen Collection and Preparation for Analysis

Routinely processed, neutral-buffered formalin-fixed, paraffin-embedded, tissues are suitable for use with this primary antibody when used with Cell Marque detection kits (see Materials, Reagents, and Equipment Needed But Not Provided section). Note: Cell Marque evaluates performance only on human tissues. The recommended tissue fixative is 10% neutral-buffered formalin. Variable results may occur as a result of prolonged fixation or special processes such as decalcification of bone marrow preparations.

Each section should be cut to the appropriate thickness (approximately 3 µm) and placed on a positively charged glass slide. Slides containing the tissue section may be baked for at least 2 hours (but not longer than 24 hours) in a 58-60°C ± 5°C oven.

## WARNINGS AND PRECAUTIONS

1. Take reasonable precautions when handling reagents. Use disposable gloves and lab coats when handling suspected carcinogens or toxic materials (example: xylene).
2. Avoid contact of reagents with eyes and mucous membranes. If reagents come in contact with sensitive areas, wash with copious amounts of water.
3. Patient specimens and all materials contacting them should be handled as biohazardous materials and disposed of with proper precautions. Never pipette by mouth.
4. Avoid microbial contamination of reagents, as this could produce incorrect results.
5. The user must validate incubation times and temperatures.
6. The prediluted, ready-to-use reagents are optimally diluted, and further dilution may result in loss of antigen staining.
7. The concentrated reagents may be diluted optimally based on validation by user. Any diluent used that is not specifically recommended herein must likewise be validated by the user for both its compatibility and effect on stability.

8. When used according to instructions, this product is not classified as a hazardous substance. The preservative in the reagent is less than 0.1% sodium azide and does not meet the OSHA (USA) criteria for hazardous substance at the stated concentration. See MSDS.
9. The user must validate any storage conditions other than those specified in the package insert.
10. Diluent may contain bovine serum albumin and supernatant may contain bovine serum. The products containing fetal bovine serum and products containing bovine serum albumin are purchased from commercial suppliers. Certificates of Origin for the animal source used in these products are on file at Cell Marque. The certificates support that the bovine sources are from countries with negligible BSE risk and state sources of bovine from USA and Canada.
11. As with any product derived from biological sources, proper handling procedures should be used.

## INSTRUCTIONS FOR USE

### Step by Step Procedure

Recommended Staining Protocols for the stated primary antibody (414R series):

#### HiDef Detection™ System:

1. Deparaffinize, rehydrate, and epitope retrieve; the preferred method is the use of Heat Induced Epitope Retrieval (HIER) techniques using Cell Marque's Trilogy™ in conjunction with a pressure cooker. The preferred method allows for simultaneous deparaffinization, rehydration, and epitope retrieval. Upon completion, rinse with 5 changes of distilled or deionized water.
2. If using HRP detection system, place slides in peroxide block for 10 minutes; rinse. If using AP detection system, omit this step.
3. Apply the antibody and incubate for 10 - 30 minutes; rinse.
4. Apply the HiDef Detection™ Amplifier Rabbit/Mouse for 10 minutes; rinse.
5. Apply the polymer detector for 10 minutes; rinse.
6. Apply ample amount of chromogen and incubate for 1 - 10 minutes; rinse.
7. Dehydrate and coverslip.

## QUALITY CONTROL PROCEDURES

### Positive Tissue Control

A positive tissue control must be run with every staining procedure performed. This tissue may contain both positive and negative staining cells or tissue components and serve as both the positive and negative control tissue. Control tissues should be fresh autopsy, biopsy or surgical specimens prepared or fixed as soon as possible in a manner identical to the test sections. Use of a tissue section fixed or processed differently from the test specimen will serve to provide control for all reagents and method steps except fixation and tissue processing.

A tissue with weak positive staining is more suitable for optimal quality control and for detecting minor levels of reagent degradation. Positive tissue control for the stated primary antibody (414R series) may include the following:

Breast carcinoma	Cytoplasmic, membranous
------------------	-------------------------

Known positive tissue controls should be utilized only for monitoring the correct performance of processed tissues and test reagents, not as an aid in determining a specific diagnosis of patient samples. If the positive tissue controls fail to demonstrate appropriate positive staining, results with the test specimens must be considered invalid.

### Negative Tissue Control

The same tissue used for the positive tissue control may be used as the negative tissue control. The variety of cell types present in most tissue sections offers internal negative control sites, but this should be verified by the user. The components that do not stain should demonstrate the absence of specific staining, and provide an indication of non-specific background staining. If specific staining occurs in the negative tissue control sites, results with the patient specimens must be considered invalid.

### Unexplained Discrepancies

Unexplained discrepancies in controls should be referred to Cell Marque Customer Service immediately. If quality control results do not meet specifications, patient results are invalid. See the Troubleshooting section of this insert. Identify and correct the problem, then repeat the entire procedure with the patient samples.

### Negative Control Reagent

A negative control reagent must be run for every specimen to aid in the interpretation of results. A negative control reagent is used in place of the primary antibody to evaluate nonspecific staining. The slide should be treated with negative control reagent, matching the host species of the primary antibody, and ideally having the same IgG concentration. The incubation period for the negative control reagent should equal the primary antibody incubation period.

## INTERPRETATION OF RESULTS

The immunostaining procedure causes a colored reaction product to precipitate at the antigen sites localized by the primary antibody. Refer to the appropriate detection system package insert for expected color reactions. A qualified pathologist experienced in immunohistochemistry procedures must evaluate positive and negative tissue controls before interpreting results.

### Positive Tissue Control

The stained positive tissue control should be examined first to ascertain that all reagents are functioning properly. The presence of an appropriately colored reaction product within the target cells is indicative of positive reactivity. Refer to the package insert of the

detection system used for expected color reactions. Depending on the incubation length and potency of the hematoxylin used, counterstaining will result in a pale to dark blue coloration of cell nuclei. Excessive or incomplete counterstaining may compromise proper interpretation of results. If the positive tissue control fails to demonstrate appropriate positive staining, any results with the test specimens are considered invalid.

#### Negative Tissue Control

The negative tissue control should be examined after the positive tissue control to verify the specific labeling of the target antigen by the primary antibody. The absence of specific staining in the negative tissue control confirms the lack of antibody cross reactivity to cells or cellular components. If specific staining occurs in the negative tissue control, results with the patient specimen are considered invalid. Nonspecific staining, if present, will have a diffuse appearance. Sporadic light staining of connective tissue may also be observed in sections from tissues that are not optimally fixed. Intact cells should be used for interpretation of staining results. Necrotic or degenerated cells show non-specific staining.

#### Patient Tissue

Patient specimens should be examined last. Positive staining intensity should be assessed within the context of any background staining of the negative reagent control. As with any immunohistochemical test, a negative result means that the antigen in question was not detected, not that the antigen is absent in the cells or tissue assayed. A panel of antibodies may aid in the identification of false negative reactions (see Summary of Expected Results section). The morphology of each tissue sample should also be examined utilizing a hematoxylin and eosin stained section when interpreting any immunohistochemical result. The patient's morphologic findings and pertinent clinical data must be interpreted by a qualified pathologist.

#### LIMITATIONS

1. This reagent is "for professional use only" as immunohistochemistry is a multiple step process that requires specialized training in the selection of the appropriate reagents, tissues, fixation, processing; preparation of the immunohistochemistry slide; and interpretation of the staining results.
2. For laboratory use only.
3. For *in vitro* diagnostic use.
4. Tissue staining is dependent on the handling and processing of the tissue prior to staining. Improper fixation, freezing, thawing, washing, drying, heating, sectioning, or contamination with other tissues or fluids may produce artifacts, antibody trapping, or false negative results. Inconsistent results may result from variations in fixation and embedding methods, as well as from inherent irregularities within the tissue.
5. Excessive or incomplete counterstaining may compromise proper interpretation of results.

6. The clinical interpretation of any positive staining, or its absence, must be evaluated within the context of clinical history, morphology, other histopathological criteria as well as other diagnostic tests. This antibody is intended to be used in a panel of antibodies if applicable. It is the responsibility of a qualified pathologist to be familiar with the antibodies, reagents, diagnostic panels, and methods used to produce the stained preparation. Staining must be performed in a certified, licensed laboratory under the supervision of a pathologist who is responsible for reviewing the stained slides and assuring the adequacy of positive and negative controls.
7. Cell Marque provides antibodies and reagents at optimal dilution for use as instructed. Any deviation from recommended test procedures may invalidate expected results. Appropriate controls must be employed and documented. Users in any circumstance must accept responsibility for interpretation of patient results.
8. Cell Marque provides primary antibodies in concentrated format so that the user may subsequently optimally dilute for use subject to the user's determination of and adherence to suitable validation techniques. Users must validate the use of any diluents other than what is recommended herein. Once the primary is validated to be suitable for use, any deviation from recommended test procedures may invalidate expected results. Appropriate controls must be employed and documented. Users in any circumstance must accept responsibility for interpretation of patient results.
9. This product is not intended for use in flow cytometry.
10. Reagents may demonstrate unexpected reactions in previously untested tissues. The possibility of unexpected reactions even in tested tissue groups cannot be completely eliminated because of biological variability of antigen expression in neoplasms, or other pathological tissues. Contact Cell Marque customer service with any suspected, documented unexpected reactions.
11. Tissues from persons infected with hepatitis B virus and containing hepatitis B surface antigen (HBsAg) may exhibit nonspecific staining with horseradish peroxidase.
12. When used in blocking steps, normal sera from the same animal source as the secondary antisera may cause false negative or false positive results because of the effect of autoantibodies or natural antibodies.
13. False positive results may be seen because of nonimmunological binding of proteins or substrate reaction products. They may also be caused by pseudoperoxidase activity (erythrocytes), endogenous peroxidase activity (cytochrome C), or endogenous biotin (example: liver, brain, breast, kidney) subject to the type of immunostaining technique used.
14. As with any immunohistochemistry test, a negative result means that the antigen was not detected, not that the antigen was absent in the cells or tissue assayed.

#### Specific Limitations

1. The prediluted antibody products are optimized as a ready-to-use product. Because of the possibility of variation in tissue fixation and processing, it may be necessary to increase or decrease the primary antibody incubation time on individual specimens.

- The antibody, in combination with detection systems and accessories, detects antigen(s) that survive routine formalin fixation, tissue processing and sectioning. Users who deviate from recommended test procedures remain, as they would in any circumstance, responsible for interpretation and validation of patient results.

### Summary of Expected Results

See the following tables of reactivity:

Normal Study			
Tissue	# Stained	Total #	Notes
Brain	0	1	
Adrenal Cortex	1	1	Medulla (-)
Ovary	1	1	
Pancreas	1	1	Islet cells (-)
Parathyroid	1	1	
Pituitary	1	1	
Testis	1	1	Leydig cells (+)
Thyroid	1	1	
Breast	1	1	
Spleen	1	1	
Tonsil	0	1	Dendritic cells (+)
Thymus	1	1	Thymus epithelia (+)
Bone Marrow	0	1	
Lung	1	1	Bronchus (+)
Heart	0	1	
Esophagus	1	1	
Stomach	0	1	
Small Intestine	1	1	
Colon	1	1	
Liver	1	1	
Salivary Gland	1	1	
Gall Bladder	1	1	
Kidney	1	1	
Bladder	1	1	
Prostate	1	1	
Uterus	1	1	
Fallopian Tube	1	1	
Ureter	1	1	
Cervix	1	1	
Skeletal Muscle	0	1	
Smooth Muscle	1	1	
Skin	1	1	
Peripheral Nerve	1	1	
Mesothelium	1	1	

Normal Study			
Fat	1	1	
Placenta	1	1	

EGFR protein expression is almost found in all normal tissue, which is indicated in the literature.

Disease Tissue Study			
Tissue	# Stained	Total #	Notes
Breast invasive ductal carcinoma	36	61	
Renal cell carcinoma	22	24	
Colorectal adenocarcinoma	19	20	
Lung adenocarcinoma	9	10	
Hepatocellular carcinoma	2	9	

This antibody performs as indicated in the literature.

### TROUBLESHOOTING

- If the positive control exhibits weaker staining than expected, other positive controls run during the same staining run should be checked to determine if it is because of the primary antibody or one of the common secondary reagents.
- If the positive control is negative, other positive controls used on the same run should be checked to determine if the underlying cause relates to the primary antibody or one of the common secondary reagents. Tissues may have been improperly collected, fixed or deparaffinized. The proper procedure should be followed for collection, storage, and fixation.
- If excessive background staining occurs, high levels of endogenous biotin may be present. A biotin blocking step should be included unless a biotin-free detection system is being used in which case any biotin present would not be a contributing factor to background staining.
- If all of the paraffin has not been removed, the deparaffinization procedure should be repeated.
- If tissue sections wash off the slide, slides should be checked to ensure that they are positively charged. Other possibilities that could have adverse affect on tissue adhesion include insufficient drying of the tissue section on the slide prior to staining or fixation in formalin that was not properly neutral-buffered. Tissue thickness may also be a contributing factor.

For corrective action, refer to the Step by Step Procedure section or contact Cell Marque customer service.



**REFERENCES**

1. Vranic, S, et al. Modern Pathology (2010) 23, 644.
2. Ch'ng, S, et al. Hum Pathol 2008;39:344.
3. Van Damme, et al. BMC Cancer 2010;10:189.
4. Tawbi, H, et al. Oncologist 2008;13:459.
5. Cascio, MJ, et al. Mod Pathol 2010; 23:574.

**DISCLAIMERS**[www.cellmarque.com](http://www.cellmarque.com)EMERGO EUROPE  
Molenstraat 15, 2513 BH, The Hague, NL.

CM Template #1.3 v1