

Anti-Influenza B Virus HA antibody, mouse monoclonal (1H12)

65-160 100 µg

Background: Hemagglutinin (HA) binds to sialic acid-containing receptors on the cell surface, bringing about the attachment of the virus particle to the cell. Plays a major role in the determination of host range restriction and virulence. Class I viral fusion protein. Responsible for penetration of the virus into the cell cytoplasm by mediating the fusion of the membrane of the endocytosed virus particle with the endosomal membrane. Low pH in endosomes induce an irreversible conformational change in HA2, releasing the fusion hydrophobic peptide. Several trimers are required to form a competent fusion pore.

Post-translational modification:

HAo consists of 584 amino acids with molecular mass of 63,275. In natural infection, inactive HA is matured into HA1 and HA2. By sequence similarity it is indicated to be palmitoylated.

Applications

- 1) Western blotting (1/500~1/1,000 dilution)
- 2) Immunofluorescent and Immunocytochemical staining (1/100~1/200 dilution)
- 3) Immunoprecipitation (1/200 dilution)
- 4) ELISA (assay dependent)

Immunogen: Human Influenza B Virus strain Nagasaki/1/87, one of the strains of B/Victoria group.

Specificity: According to Ref.1 during epidemic in Osaka 1996-97, 1H12 antibody reacted with HA protein of most of Influenza B virus isolates belonging to Victoria group tested (70/73 strains) and none of clinical 27 isolates belonging to Yamagata group as examined by PAP staining. However later test with IF staining it also reacts with some of Yamagata group strains like B/Florida/4/2006. However, note that HA changes during passages and may change reactivity to this antibody.

No cross reactivity with any strains of influenza A virus.

Isotype: mouse IgG1, kappa

Purity: Produced in serum-free medium and purified by proprietary chromatography procedure under mild conditions.

Form: 1 mg/ml in PBS, 50% glycerol, filter sterilized

Storage: Shipped at 4°C or -20°C, and upon arrival, store at -20°C.

Data Link: UniProtKB [P03460](#) Influenza B/Lee/1940 HA protein.



COSMO BIO CO., LTD.

TOKYO JAPAN
TEL : (81)3-5632-9617
FAX : (81)3-5632-9618
e-mail : export@cosmobio.co.jp



www.cosmobio.com

References: This antibody was described and used in the following reference.

1. Nakagawa N. et al. Rapid detection and identification of two lineages of influenza B strains with monoclonal antibodies. [J Virol Methods](#). 1999;79:113-2 ICC, IP

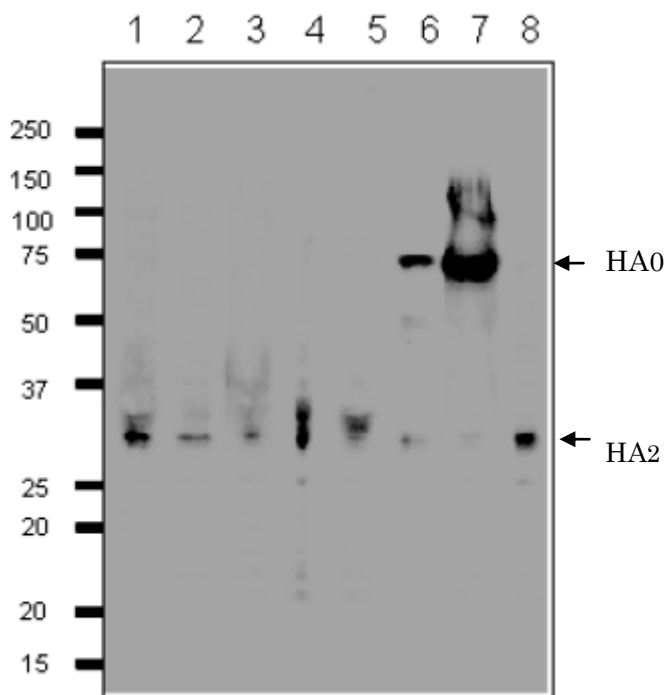


Fig. 1 Detection of HA protein in the crude extracts of MDCK cells infected with various Influenza B virus strains by western blotting using 1H12 monoclonal antibody.

1. B/Mie/1/1993
2. B/Johannes Burg/5/1999
3. B/Florida/4/2006
4. B/Lee/1940
5. B/Florida/4/2006
6. Shandong/7/97
7. B/Malaysia/2506/2004
8. B/Massachusetts/2/2012

First antibody was used at 1/500 dilution and as 2nd antibody, HRP-conjugated goat anti-mouse IgG antibody was used at 1/10,000 dilution. Positions of marker proteins are indicated in kDa on the left. **Clone 1H12 recognizes an epitope on HA2 region..**

Mock

B/Malaysia/2506/2004

B/Florida/4/2006



COSMO BIO CO., LTD.

TOKYO JAPAN

TEL : (81)3-5632-9617

FAX : (81)3-5632-9618

e-mail : export@cosmobio.co.jp



www.cosmobio.com

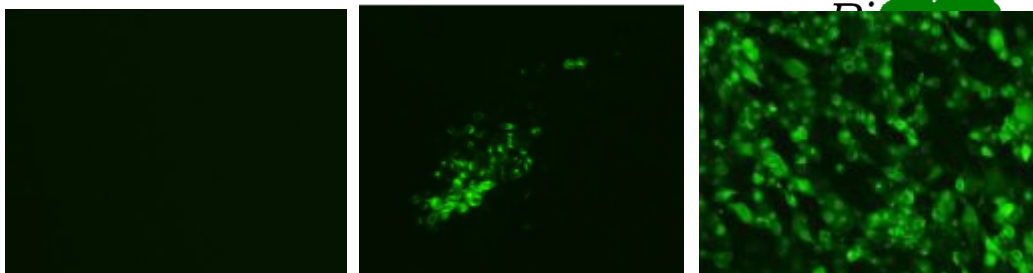


Fig.2 Immunofluorescence assay of MDCK (canine kidney) cells infected with Influenza B virus, using anti-Influenza B virus HA antibody (clone 1H12).

Anti-Influenza B Virus HA antibody (clone 1H12) efficiently detected the viruses in the infected MDCK cells with B/Malaysia/2506/2004 and B/Florida/4/2006 virus strains. The cells were fixed with 4% paraformaldehyde in phosphate-buffered saline (PBS) and permeabilized with 0.1% Triton X-100 in PBS. The bound antibody was visualized by a further reaction with an Alexa Fluor 488-conjugated secondary antibody (green). Image on the left is a negative control, mock-infected MDCK cells.



COSMO BIO CO., LTD.

TOKYO JAPAN
TEL : (81)3-5632-9617
FAX : (81)3-5632-9618
e-mail : export@cosmobio.co.jp



www.cosmobio.com