

POLYCLONAL ANTIBODY

For research use only. Not for clinical diagnosis.

Catalog No. CNP- IP-238

Anti TALDO1

Product type Primary antibodies

Immunogen TALDO1 (NP_006746.1, 1 a.a. ~ 337 a.a) full-length human protein.

Form Liquid
Protein A purified

Volume 100ul

Specificity Human


Host Rabbit

Immunogen Sequence MSSSPVKRQRMESALDQLKQFTTVVADTGDFHAIDEYKPQDATTNPSLILAAAQMPAYQE
LVEEAIAYGRKLGGSQEDQIKNAIDKLFVLFGAIEILKKIPGRVSTEVDARLSFDKDAMVA
RARRLIELYKEAGISKDRILIKLSSTWEGIQAGKELEEQHGHCNMTLLFSFAQAVACAE
AGVTLISPFVGRILDWHVANTDKKSYEPLDPGVKSVTKIYNYYKKFSYKTIVMGASFRN
TGEIKALAGCDFLTISPKLLGELLQDNAKLVPVLSAKAAQASDLEKIHLDKESFRWLHNE
DQMAVEKLSDGIRKFAADAVKLERMLTERMFNAENGK

Gene Alias TAL, TAL-H, TALDOR, TALH

Application notes **WB**
Antibody dilution: 1:1000. Western Blot analysis of TALDO1 expression in HepG2 (Predicted Mw: 37.5400 kDa).

Western Blotting (Transfected lysate)

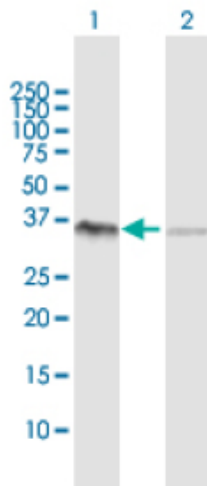


Western Blot analysis of TALDO1 expression in transfected 293T cell line by TALDO1 rabbit polyclonal antibody.

Lane 1: TALDO1 transfected lysate(37.50 KDa).

Lane 2: Non-transfected lysate.

Antibody dilution: 1:1000. Western Blot analysis of TALDO1 expression in human pancreas (Predicted Mw: 37.5400 kDa).

Western Blotting (Tissue lysate)

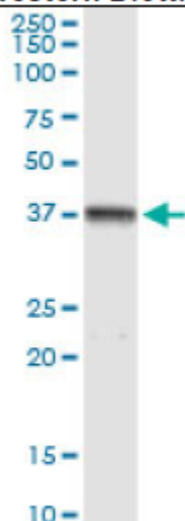


TALDO1 rabbit polyclonal antibody. Western Blot analysis of TALDO1 expression in human pancreas.

IP

Antibody dilution 1:15 ~ 1:60 for immunoprecipitation of TALDO1 transfected lysate.

Immunoprecipitation

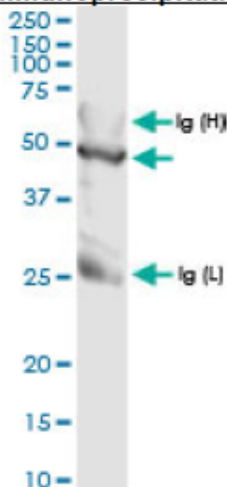

Immunoprecipitation of TALDO1 transfected lysate using TALDO1 rabbit polyclonal antibody and Protein A Magnetic Bead, and immunoblotted with TALDO1 rabbit polyclonal antibody.

Storage

For continuous use, store at 2-8°C for one-two days. For extended storage, store in -20°C freezer. Working dilution samples should be discarded if not used within 12 hours.
For research use only. Not for clinical diagnosis.



COSMO BIO Co., LTD.

Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

URL: <http://www.cosmobio.co.jp>

e-mail: export@cosmobio.co.jp

[Outside Japan] Phone : +81-3-5632-9617

[国内連絡先] Phone : +81-3-5632-9610

FAX : +81-3-5632-9618

FAX : +81-3-5632-9619