

POLYCLONAL ANTIBODY

For research use only. Not for clinical diagnosis.

Catalog No. CNP-STOML2-415

Anti STOML2

Product type Primary antibodies

Immunogen STOML2 (NP_038470.1, 1 a.a. ~ 356 a.a) full-length human protein.

Form Liquid
Protein A purified

Volume 100ul

Specificity Human


Host Rabbit

Immunogen Sequence MLARAARGTGALLLRGSLASGRAPRRASSGLPRNTVVLFPQQEAWVVERMGRFHRILE
PGLNILIPVLDRIYVQSLKEIVINVEQSAVTLDNVTLQIDGVLYLRIMDPYKASYGVE
DPEYAVTQLAQTTMRSELGKLSLDKVFREERESLNASIVDAINQAADCWGI RCLRYEIKDI
HVPPRVKESMQMQVEAERRKRATVLESEGTRESAINVAEGKKQAQILASEAEKAEQINQA
AGEASAVLAKAKAKAEAIRILAAALTQHNGDAAASLTVAEQYVSAFSKLAKDSNTILLPS
NPGDVTSMVAQAMGVYGALTKAPVPGTPDSLSSGSSRDVQGT DASLDEELDRVKMS

Gene Alias HSPC108, SLP-2

Application notes **WB**
Antibody dilution: 1:1000 ~ 1:2000. Specific band of ~38.50 kDa using STOML2 transfected 293T lysate.

Western Blotting (Transfected lysate)

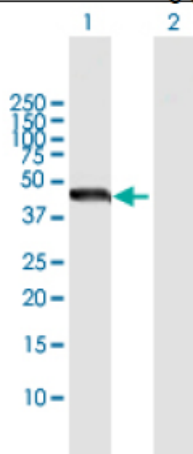


Western Blot analysis of STOML2 expression in transfected 293T cell line by STOML2 rabbit polyclonal antibody.

Lane 1: STOML2 transfected lysate(38.50 KDa).

Lane 2: Non-transfected lysate.

Antibody dilution: 1:1000. Western Blot analysis of STOML2 expression in HepG2 (Predicted Mw: 38.5340 kDa).

Western Blotting (Cell lysate)

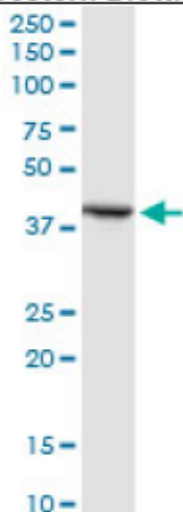


STOML2 rabbit polyclonal antibody. Western Blot analysis of STOML2 expression in HepG2.

IP

Antibody dilution 1:15 ~ 1:60 for immunoprecipitation of STOML2 transfected lysate.

Immunoprecipitation

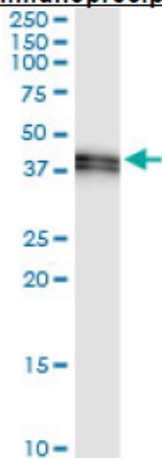

Immunoprecipitation of STOML2 transfected lysate using STOML2 rabbit polyclonal antibody and Protein A Magnetic Bead, and immunoblotted with STOML2 purified mouse polyclonal antibody.

Storage

For continuous use, store at 2-8°C for one-two days. For extended storage, store in -20°C freezer. Working dilution samples should be discarded if not used within 12 hours.

For research use only. Not for clinical diagnosis.



COSMO BIO Co., LTD.

Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

URL: <http://www.cosmobio.co.jp>

e-mail: export@cosmobio.co.jp

[Outside Japan] Phone : +81-3-5632-9617

[国内連絡先] Phone : +81-3-5632-9610

FAX : +81-3-5632-9618

FAX : +81-3-5632-9619