

MONOCLONAL ANTIBODY

For research use only, Not for diagnostic use.

Catalog No. CAC-SHI-EXO-M01-EX

Anti CD9 [Clone : 12A12]

BACKGROUND

CD9 is a cell surface glycoprotein which belongs to the tetraspanin superfamily. CD9 is known to complex with integrins and other transmembrane 4 superfamily proteins. It is found on the surface of exosomes. It can modulate cell adhesion and migration and also trigger platelet activation and aggregation.

Product type Primary antibody

Immunogen Human CD9(a.a 113-195)

Rased in Mouse Myeloma -

Clone number 12A12 Isotype IgG2b

Source Culture supernatant

Purification Affinity purified by Protein A

Buffer Phosphate buffered saline containing 0.1% NaN₃ as a preservative

Concentration1 mg / mLVolume100 ulLabelUnlabeledSpecificityHuman CD9

Cross reactivity Human. Other species is not tested.

Storage Store below - 40°C

Aliquot to avoid cycles of freeze/thaw..

Application notes
Recommended
dilutions

• Western blotting: 1/500 - 1/2,000 (0.5 - 2ug/mL, non-reductive condition)

Immunoprecipitation: 1/200 – 1/1,000 (1- 5ug/mL)

* Experiment example : Serum Exosome

ELISA: 1/10,000 – 1/200,000 (5 ng/mL - 0.1ug/mL)

Other applications have not been tested.

Optimal dilutions/concentrations should be determined by the end user.

References

ANTIBODY CHARACTERIZATION

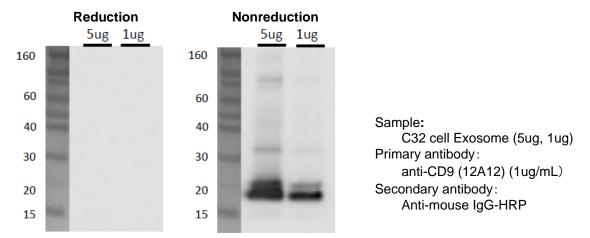


Fig.1 Western blot analysis of CD9 in C32 cell Exosome with anti-CD9 (12A12) antibody



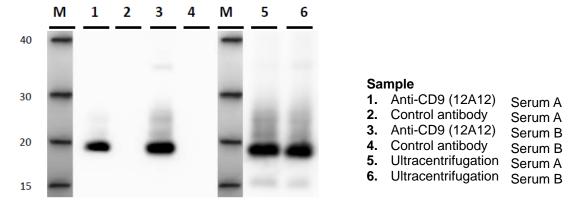


Fig.2 Immunoprecipitation in Human serum exosome

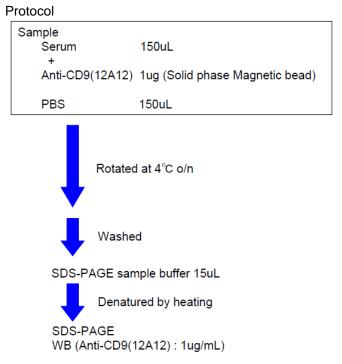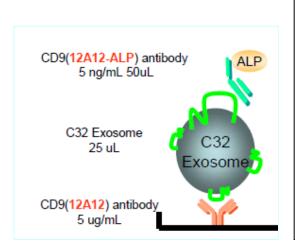

ANTIBODY CHARACTERIZATION



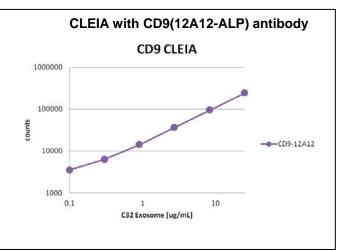


Fig.3 ELISA data of C32 cell exosome with anti-CD9 (12A12-ALP) antibody (5 ng/mL)

RELATED PRODUCT:

Product Name	Application	Quantity	Cat#
Anti CD9 (12A12) Monoclonal Antibody	WB / IP / ELISA	100 ug	CAC-SHI-EXO-M01-EX
Anti CD63 (8A12) Monoclonal Antibody	WB / IP	100 ug	CAC-SHI-EXO-M02-EX
Anti CD81 (12C4) Monoclonal Antibody	WB / IP	100 ug	CAC-SHI-EXO-M03-EX

For research use only, Not for diagnostic use.





TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

URL: http://www.cosmobio.co.jp e-mail: export@cosmobio.co.jp

[Outside Japan] Phone: +81-3-5632-9617 [国内連絡先] Phone: +81-3-5632-9610 FAX: +81-3-5632-9618 FAX: +81-3-5632-9619