




Code No.KAL-KR099

For research use only

Anti Human 17 β HSD 8 Polyclonal Antibody

17 β HSD 8 : 17 β hydroxysteroid dehydrogenase type 8

Package Size	25 μ g (100 μ L / vial)
Format	Rabbit polyclonal antibody 0.25mg/mL
Buffer	PBS [containing 2% Block Ace as a stabilizer, 0.1% Proclin as a bacteriostat]
Storage	Store below -20°C Once thawed, store at 4°C. Repeated freeze-thaw cycles should be avoided.
Purification method	This antibody was purified from rabbit serum immunized with synthesized partial peptide of human 17 β HSD8 by peptide affinity chromatography.
Working dilution	For Immunohistochemistry ; 5 μ g/mL

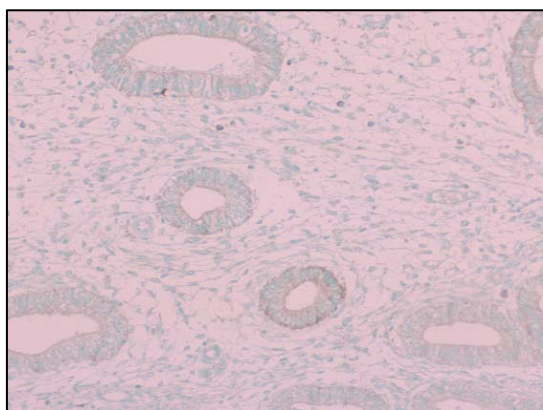

Immunohistochemistry

Sample : Human endometrium

Preparation of antibodies and instruction :

Dr. Okamura H. at Department of Reproductive Medicine and Surgery, Faculty of Medical and Pharmaceutical Sciences
Kumamoto University

Distributor



COSMO BIO CO., LTD.
Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

<http://www.cosmobio.co.jp>

Phone : +81-3-5632-9617

e-mail : export@cosmobio.co.jp

FAX : +81-3-5632-9618



Code No.KAL-KR099

Anti Human 17 β HSD 8 Polyclonal Antibody

17 β HSD 8 : 17 β hydroxysteroid dehydrogenase type 8

【Reference】

1. Ramirez S. et al. : Mol Cell Endocrinol. 1998 Aug 25;143(1-2):9-22
2. Woo D. et al. : Mol Cell Endocrinol. 2001 May 15;176(1-2):155-62

Distributor



COSMO BIO CO., LTD.
Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

<http://www.cosmobio.co.jp>

Phone : +81-3-5632-9617


e-mail : export@cosmobio.co.jp

FAX : +81-3-5632-9618

抗 ヒト 17 β HSD 8 ポリクローナル抗体

17 β HSD 8 : 17 β hydroxysteroid dehydrogenase type 8

容量	25 μ g (100 μ L/vial)
形状	ウサギポリクローナル抗体 0.25mg/mL、凍結品
バッファー	PBS [2%ブロッカーエース(安定化蛋白)、0.1%proclin 含有]
保管方法	-20℃以下 抗体を低濃度にて冷蔵保管されますと、失活する恐れがあります。 融解後は4℃で保存し、お早めにご使用下さい。 また凍結融解を繰り返すことは避けて下さい。
製造方法	ヒト 17 β HSD 8 の部分ペプチドを免疫して得られたウサギの抗血清より、ペプチドアフィニティーカラムにて精製。
使用濃度	免疫組織染色: 5 μ g/mL

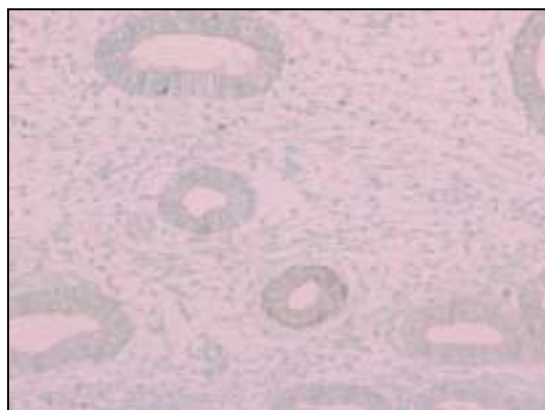

免疫組織染色

Sample : ヒト子宮内膜
(増殖期子宮内膜に陽性)

提供 : 熊本大学大学院 医学薬学研究部
先端生命医療科学部門 成育医学講座
岡村 均 先生
永吉 裕三子 先生



人と科学のステキな未来へ

コスモ・バイオ株式会社

〒135-0016 東京都江東区東陽 2-2-20 東陽駅前ビル
URL : <http://www.cosmobio.co.jp/>

● 営業部 (お問い合わせ)

TEL : (03) 5632-9610 FAX : (03) 5632-9619
TEL : (03) 5632-9620

抗 ヒト 17β HSD 8 ポリクローナル抗体

17β HSD 8 : 17β hydroxysteroid dehydrogenase type 8

【参考文献】

1. Ramirez S. et al. : Mol Cell Endocrinol. 1998 Aug 25;143(1-2):9-22
2. Woo D. et al. : Mol Cell Endocrinol. 2001 May 15;176(1-2):155-62



人と科学のステキな未来へ

コスモ・バイオ株式会社

〒135-0016 東京都江東区東陽 2-2-20 東陽駅前ビル
URL : <http://www.cosmobio.co.jp/>

● 営業部（お問い合わせ）

TEL : (03) 5632-9610 FAX : (03) 5632-9619

TEL : (03) 5632-9620