



## Anti PACAP 27 (Human) Serum

Cat. No. YII-Y050-EX      Lot No. 080271013

**Description:** This antiserum was raised in a rabbit by immunization with a carrier free synthetic PACAP 27 (human) peptide. The product vial contains 50 µL of the titled antiserum obtained by lyophilizing its 0.001 M phosphate buffer (pH 7.0, 0.5mL) solution. It can be used for immunoassay, immunohistochemistry or any other immunoreaction with PACAP27 (human).

**Immunogen:** Synthetic PACAP 27 (human) peptide, carrier free **Host:** Rabbit

### **Amino Acid Sequence of PACAP27 (human)<sup>1)</sup>:**

HSDGIFTDSY SRYRKQMAVK KYLA AVL-NH<sub>2</sub>

**Product Form:** Lyophilized unpurified serum **Size:** 50 µL

**Reconstitution:** Reconstitute the product with 0.5mL of 0.01M PBS (pH 7.0) to make a 10 fold diluted stock solution. If it is stored in a refrigerator, add moderate antiseptic to the solution (e.g. NaN<sub>3</sub> 0.1%).

**Storage:** The product will be stable for over one year if it be stored at -20°C to -80°C until opened. Upon reconstitution, the antiserum solution must be stored at 2°C to 8°C and used within one month. Repeated freezing-thawing should be avoided.

**Suggested Working Dilution Range:** 1:5,000-40,000 (final dilution ~1:280,000) for radioimmunoassay;

1: 1,000-4,000 for immunohistochemistry (frozen or paraffin sections). Optimal dilution should be determined by each laboratory for each application.

**Specificity** (based on radioimmunoassay): PACAP27 (human) 100%, PACAP38 (human) 0%, VIP (porcine) 0%, PHI (rat) 0%, PHI (human) 0%, secretin (human) 0%, helodermin 0%

**Positive Control** (immunohistochemistry): Rat colon

**Species Tested:** Rat

### **REFERENCES:**

1) A. Miyata, A. Arimura et al., Isolation of a neuropeptide corresponding to the N-terminal 27 residues of the pituitary adenylate cyclase activating polypeptide with 38 residues (PACAP38). *Biochemistry Biophysics Research Communications* 170:643-648, 1990

### **FOR RESEARCH LABORATORY USE ONLY**

DO NOT USE ORGANIC SOLVENTS FOR DISSOLVING ANTISERUM


#### **Distributor**



COSMO BIO CO., LTD.  
Inspiration for Life Science

TOYO EKIMAE BLDG. 2-20, TOYO 2CHOME  
KOTO-KU, TOKYO 135-0016, JAPAN  
TEL : +81-3-5632-9617  
FAX : +81-3-5632-9618  
URL : <http://www.cosmobio.co.jp/>  
e-mail : [export@cosmobio.co.jp](mailto:export@cosmobio.co.jp)

#### **Manufacturer**

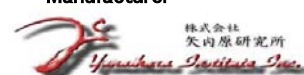

2480-1 Awakura, Fujinomiya-shi,  
Shizuoka 418-0011, JAPAN  
Phone: +81-544-22-2771  
Fax: +81-544-22-2770  
URL: [www.yanaihara.co.jp](http://www.yanaihara.co.jp)  
E-mail: [ask@yanaihara.co.jp](mailto:ask@yanaihara.co.jp)