

POLYCLONAL ANTIBODY

For research use only. Not for clinical diagnosis.

Catalog No. CNP-FANCL-075


Anti FANCL

Product type	Primary antibodies
Immunogen	FANCL (NP_060532.2, 1 a.a. ~ 375 a.a) full-length human protein.
Form	Liquid
	Protein A purified
Volume	100ul
Specificity	Human
Host	Rabbit

Immunogen Sequence	MAVTEASLLRQCPLLLPQNRSKTVYEGFISAQGRDFHLRIVLPEDLQLKNARLLCSWQLR TILSGYHRIVQQRMQHSPDLMSEFMELKMLLEVALKNRQELYALPPPQFYSSLIEEIGT LGWDKLVYADTCFSTIKLKAEDASGREHLITLKLKAKYPAESPDYFVDFPVPFCASWTPQ SSLISISQFLAAIESLKAFWDVMDIDEKWTWLEPEKPPRSATARRIALGNNVSINIEV DPRHPTMLPECFGLGADHVVKPLGIKLSRNIHLWDPENSVLQNLKDVLEIDFPARAILEK SDFTMDCGICYAYQLDGTIPDQVCDNSQCGQPFHQICLYEWLRGLLTSRQSFNIIFGECF YCSKPITLKMSGRKH
---------------------------	--

Gene Alias	FAAP43, FLJ10335, PHF9, POG
Application notes	WB Antibody dilution: 1:1000 ~ 1:2000. Specific band of ~42.90 kDa using FANCL transfected 293T lysate.

Western Blotting (Transfected lysate)

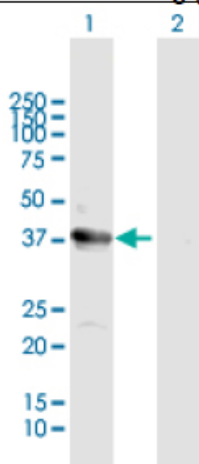


Western Blot analysis of FANCL expression in transfected 293T cell line by FANCL rabbit polyclonal antibody.

Lane 1: FANCL transfected lysate(42.90 kDa).

Lane 2: Non-transfected lysate.

Antibody dilution: 1:1000. Western Blot analysis of FANCL expression in human spleen (Predicted Mw: 42.9150 kDa).

Western Blotting (Tissue lysate)

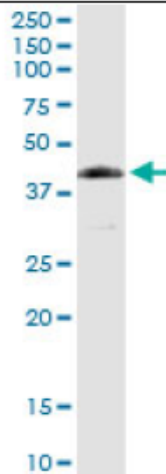


FANCL rabbit polyclonal antibody. Western Blot analysis of FANCL expression in human spleen.

IP

Antibody dilution 1:15 ~ 1:60 for immunoprecipitation of FANCL transfected lysate.

Immunoprecipitation

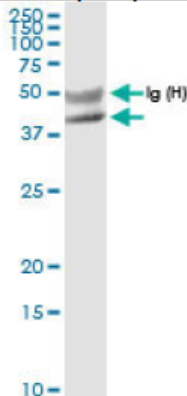

Immunoprecipitation of FANCL transfected lysate using FANCL rabbit polyclonal antibody and Protein A Magnetic Bead, and immunoblotted with FANCL rabbit polyclonal antibody.

I

Storage

For continuous use, store at 2-8°C for one-two days. For extended storage, store in -20°C freezer. Working dilution samples should be discarded if not used within 12 hours.

For research use only. Not for clinical diagnosis.



COSMO BIO Co., LTD.

Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

URL: <http://www.cosmobio.co.jp>

e-mail: export@cosmobio.co.jp

[Outside Japan] Phone : +81-3-5632-9617

[国内連絡先] Phone : +81-3-5632-9610

FAX : +81-3-5632-9618

FAX : +81-3-5632-9619