

Catalog No.NU-01-PLC**Anti plectin monoclonal antibody(C-terminal)**

Product type	Primary Antibodies
Immunogen	An expressed recombinant His-tagged fusion protein of human Plectin (U53204, 2,930 - 3,153aa).
Clone number	PC742
Isotype	IgG1
Host	Mouse
Formulation	Hybridoma supernatant with 0.02% NaN ₃ as a preservative.
Volume	500ul
Label	Unlabeled
Specificity	C-terminal region of plectin
Cross reactivity	Human, mouse, rabbit, pig, bovine
Storage	Store at -20°C or -70°C in small aliquots for prolonged storage. Repeated freeze-thaw cycles can damage the immunoreactivity of an antibody.

Application notes

Recommended use	WB, IF, IP Not tested yet in other applications.
Recommended dilutions	Western Blot: 1:50–1:200 for detection of about 500 kDa polypeptide in keratinocyte cell lysates Immunohistochemistry: 1:50-1:200 for staining of acetone-fixed cryostat tissue sections. Optimal dilutions must be determined by end user.

References	Hirako Y, Yonemoto Y, Yamauchi T, Nishizawa Y, Kawamoto Y, Owaribe K. Isolation of a hemidesmosome-rich fraction from a human squamous cell carcinoma cell line. Exp. Cell Res., 324:172-182, 2014
-------------------	--

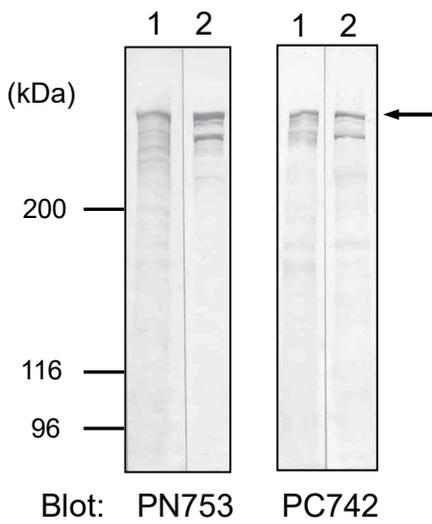


Fig.1 Western blot analysis

Whole cell lysaes prepared from DJM-1 cells (lane 1) and HeLa cells (lane 2) were immunoblotted with PN753 or PC742 at 1:200 dilution. Plectin antibodies detected approximate 500 kDa bands in these cell lysates (arrow). Smaller polypeptide found in lane 2 may be a degraded product or alternatively spliced rod-less isoform of plectin. Polypeptiedes were separated by SDS-PAGE (5% separating gel).



Fig.2 Location of the epitopes for the plectin antibodies PN753 and PC742 clones were obtained by immunizing mice with the NH2- (173-595aa) or the COOH-terminal (2,930-3,153aa) regions of human plectin (4,574aa), respectively. Gray box represents a predicted coiled-coil region (1,300-2,600aa).

Fig.3 Immunofluorescence microscopy of human skin Human skin sections were stained with PC742 (1:100 dilution) or PN753 (1:100 dilution). Arrows indicate dermal-epidermal junctions. PN753 stains epidermal cells in addition to hemidesmosomes at the dermal-epidermal junction. Sections were fixed with -20°C acetone for 10 min. Bar: 50um.

For research use only, Not for diagnostic use.



COSMO BIO Co., LTD.

【JAPAN】
TOYO EKIMAE BLDG. 2-20, TOYO 2-CHOME,
KOTO-KU. TOKYO 135-0016, JAPAN
Phone: +81-3-5632-9610
FAX: +81-3-5632-9619
URL: <https://www.cosmobio.co.jp/>



COSMO BIO USA

【Outside Japan】
2792 Loker Ave West, Suite 101
Carlsbad, CA 92010, USA
email: info@cosmobioussa.com
Phone/FAX: (+1) 760-431-4600
URL: www.cosmobioussa.com